

Classic	Van Gogh	Rembrandt
A Swivel radius*	645-715mm	665-735mm
B Minimum wall to seat back	100mm	120mm
C Seat depth	370mm	370mm
D Minimum folded width to footplate	404mm	424mm
E Minimum open width to edge of footplate	624mm	644mm
F Armrest width – external	555mm	555mm
G Armrest width – internal	439mm	439mm
H Seat back height	430mm	430mm
I Footplate to seat height	452mm	452mm
J Minimum footplate height	133mm	133mm
Bottom parking space required	650mm	650mm
Minimum track intrusion into staircase	170mm	190mm
Maximum weight limit	125kg	125kg
Speed	0.11m/sec (approx)	
Power supply	24v DC battery	
Charger	100-240v (50-60Hz)	
Maximum staircase angle	60°	

*Dimensions vary depending on hand of seat



User support: 08700 118282

Handicare, 82 First Avenue, Pensnett Estate, Kingswinford, West Midlands DY6 7FJ, United Kingdom

Textphone: (180015) 020 7490 8913 ext 05 Voicephone: (18002) 020 7490 8913 ext 05

Global Dealer Support: +44 (0) 844 225 3121 enquiries@handicare.co.uk www.handicare.co.uk

Please be aware that actual colours may vary slightly from those shown. All information is provided in good faith but does not constitute a guarantee. Handicare reserves the right to alter designs and specifications without prior notice.



handicare



Classic seat

Single tube curved stairlift system



handicare

The benefits of the Classic seat and stairlift system from Handicare

If you are thinking about installing a stairlift, and your stairs turn corners or are spiral, Handicare's Classic seat and single tube stairlift system offers you a safe and reliable solution.

Handicare stairlifts have helped thousands of people overcome the problems posed by going up and down their stairs at home.

Handicare's single tube stairlift system suits most staircases with bends or corners and when installed with the Classic seat it offers you all the basic functions that you would expect plus the option to upgrade to a powered swivel and/or have the fold up footplate linked to the seat so that folding the seat will automatically fold the footplate up too.


With a choice of three leather look upholstery colours and three standard track colours the Classic seat single tube stairlift lift system will compliment most colour schemes.

The compact unit folds up leaving the staircase free for other users and the bespoke track is designed to fit your individual staircase and maximise clear space on the stairway. It can be supplied with a hinged track if you need to avoid obstructing doorways close to the bottom of your stairs.



Fitted to the stairs, not the wall, it is unlikely you will need to redecorate following the installation of your Handicare stairlift.

Your Handicare stairlift will also come with two remote controls which means that, once folded, you can easily send it to the upstairs charge point freeing up space in your hallway.

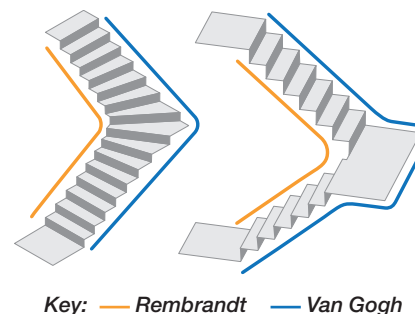

Which track?

Every staircase that is curved or turns a corner will have a shorter and a longer side. You may hear stairlift tracks being referred to as inside and outside or internal and external depending on which side to the staircase they sit on.

Handicare's single tube system has two names: the 'Rembrandt' travels along the shorter side of your staircase (the internal) and the 'Van Gogh' travels the longer side.

Your stairlift supplier will advise on which is the best solution for your staircase.

Safety as standard




Your Handicare stairlift will come with safety sensors that are designed to stop your stairlift if it meets any obstruction as you travel up the stairs and will come with a retractable seatbelt as standard.

Easy to use

The operating joystick is easy to use, simply hold it in the direction you wish to travel, up or down, and your Handicare stairlift will move smoothly and gently on its way.



All Handicare stairlifts have a battery backup which means that your stairlift will continue to work in the event of power cut so you do not need to worry about being stranded up or downstairs.

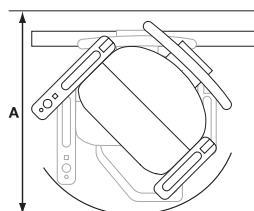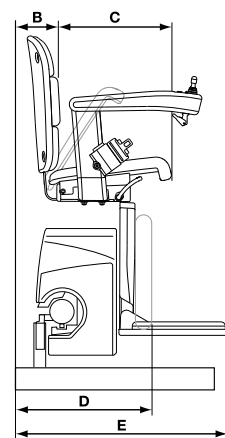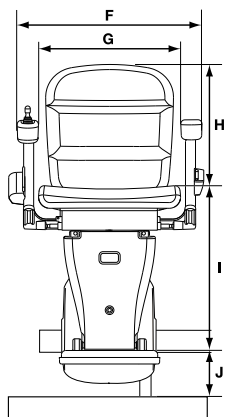

Elegance	Van Gogh	Rembrandt
A Swivel radius*	685-730mm	705-750mm
B Minimum wall to seat back	100mm	120mm
C Seat depth	390mm	390mm
D Minimum folded width to footplate	410mm	430mm
E Minimum open width to edge of footplate	620mm	640mm
F Armrest width – external	570-590mm	570-590mm
G Armrest width – internal	440-460mm	440-460mm
H Seat back height	420mm	420mm
I Footplate to seat height	495mm	495mm
J Minimum footplate height	133mm	133mm
Bottom parking space required	650mm	650mm
Minimum track intrusion into staircase	170mm	190mm
Maximum weight limit	125kg	125kg
Speed	0.11m/sec (approx)	
Power supply	24v DC battery	
Charger	100-240v (50-60Hz)	
Maximum staircase angle	60°	60°

*Dimensions vary depending on hand of seat



User support: 08700 118282

Handicare, 82 First Avenue, Pensnett Estate, Kingswinford, West Midlands DY6 7FJ, United Kingdom

Textphone: (180015) 020 7490 8913 ext 05 Voicephone: (18002) 020 7490 8913 ext 05

Global Dealer Support: +44 (0) 844 225 3121 enquiries@handicare.co.uk www.handicare.co.uk

Please be aware that actual colours may vary slightly from those shown. All information is provided in good faith but does not constitute a guarantee. Handicare reserves the right to alter designs and specifications without prior notice.



handicare



Elegance seat

Single tube curved stairlift system



handicare

The benefits of the Elegance seat and stairlift system from Handicare

If you are thinking about installing a stairlift, and your stairs turn corners or are spiral, Handicare's Elegance seat and single tube stairlift system offers you a safe and reliable solution.

Handicare stairlifts have helped thousands of people overcome the problems posed by going up and down their stairs at home.

Handicare's single tube stairlift system suits most staircases with bends or corners and when installed with the Elegance seat it offers you all the basic functions that you would expect plus the option to upgrade to a powered swivel and/or have the fold up footplate linked to the seat so that folding the seat will automatically fold the footplate up too.

With a choice of three leather look padded upholstery colours and three standard track colours the Elegance seat single tube stairlift lift system will compliment most colour schemes.




All Handicare stairlifts have a battery backup which means that your stairlift will continue to work in the event of power cut so you do not need to worry about being stranded up or downstairs.

The compact unit folds up leaving the staircase free for other users and the bespoke track is designed to fit your individual staircase and maximise clear space on the stairway. It can be supplied with a hinged track if you need to avoid obstructing doorways close to the bottom of your stairs.



Fitted to the stairs, not the wall, it is unlikely you will need to redecorate following the installation of your Handicare stairlift.

Your Handicare stairlift will also come with two remote controls which means that, once folded, you can easily send it to the upstairs charge point freeing up space in your hallway.

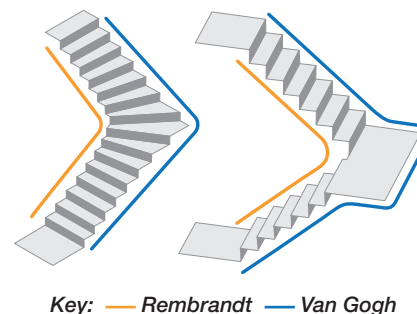

Key: — Rembrandt — Van Gogh



Which track?

Every staircase that is curved or turns a corner will have a shorter and a longer side. You may hear stairlift tracks being referred to as inside and outside or internal and external depending on which side to the staircase they sit on.

Handicare's single tube system has two names: the 'Rembrandt' travels along the shorter side of your staircase (the internal) and the 'Van Gogh' travels the longer side.




Your stairlift supplier will advise on which is the best solution for your staircase.

Safety as standard

Your Handicare stairlift will come with safety sensors that are designed to stop your stairlift if it meets any obstruction as you travel up the stairs and will come with a retractable seatbelt as standard.


Easy to use

The operating joystick is easy to use, simply hold it in the direction you wish to travel, up or down, and your Handicare stairlift will move smoothly and gently on its way.

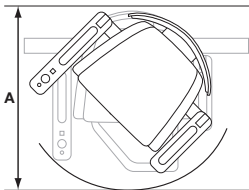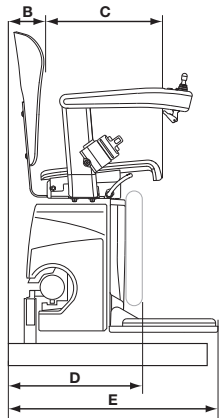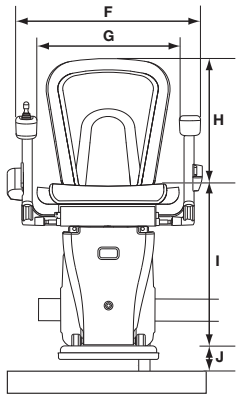

Alliance	Van Gogh	Rembrandt
A Swivel radius*	685-730mm	705-750mm
B Minimum wall to seat back	100mm	120mm
C Seat depth	390mm	
D Minimum folded width to footplate	410mm	430mm
E Minimum open width to edge of footplate	620mm	640mm
F Armrest width – external	550-590mm	
G Armrest width – internal	420-460mm	
H Seat back height	420mm	
I Footplate to seat height	485mm	
J Minimum footplate height In case of drop nose	130mm	90mm
Bottom parking space required	650mm	
Minimum track intrusion into staircase	170mm	190mm
Maximum weight limit	125kg	
Speed	0.11m/sec (approx)	
Power supply	24v DC battery	
Charger	100-240v (50-60Hz)	
Maximum staircase angle	60°	

*Dimensions vary depending on hand of seat

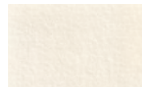

CREAM

User support: 08700 118282

Handicare, 82 First Avenue, Pensnett Estate, Kingswinford, West Midlands DY6 7FJ, United Kingdom
Textphone: (180015) 020 7490 8913 ext 05 Voicephone: (18002) 020 7490 8913 ext 05
Global Dealer Support: +44 (0) 844 225 3121 enquiries@handicare.co.uk www.handicare.co.uk

Please be aware that actual colours may vary slightly from those shown. All information is provided in good faith but does not constitute a guarantee. Handicare reserves the right to alter designs and specifications without prior notice.



Alliance seat

Single tube curved stairlift system



handicare

The benefits of the Alliance seat and stairlift system from Handicare

If you are thinking about installing a stairlift, and your stairs turn corners or are spiral, Handicare's Alliance seat and single tube stairlift system offers you a safe and reliable versatile solution.

Handicare stairlifts have helped thousands of people overcome the problems posed by going up and down their stairs at home.

Handicare's single tube stairlift system suits most staircases with bends or corners. When installed with the stylish Alliance seat it offers you all the standard functions that you would expect plus the choice to upgrade to a powered swivel and/or have the fold up footplate linked to the seat so that folding the seat will automatically fold the footplate too. This means you do not need to bend to fold and unfold your footplate.

With Alliance you can also select an offset footplate and possibility of riding while the seat is in a swivelled position. This flexibility makes for a versatile system that can accommodate compact curved staircases.

With its cream leather look upholstery colour, dark wood surround and three standard track colours the Alliance seat single tube stairlift lift system will compliment most colour schemes.

All Handicare stairlifts have a battery backup which means that your stairlift will continue to work in the event of power cut so you do not need to worry about being stranded up or downstairs.


The compact unit folds up leaving the staircase free for other users and the bespoke track is designed to fit your individual staircase and maximise clear space on the stairway.

It can be supplied with a hinged track if you need to avoid obstructing doorways close to the bottom of your stairs.



Fitted to the stairs, not the wall, it is unlikely you will need to redecorate following the installation of your Handicare stairlift.

Your Handicare stairlift will also come with two remote controls which means that, once folded, you can easily send it to the upstairs charge point freeing up space in your hallway.

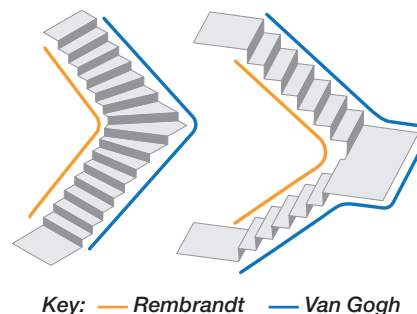

Key: — Rembrandt — Van Gogh



Which track?

Every staircase that is curved or turns a corner will have a shorter and a longer side. You may hear stairlift tracks being referred to as inside and outside or internal and external depending on which side to the staircase they sit on.

Handicare's single tube system has two names: the 'Rembrandt' travels along the shorter side of your staircase (the internal) and the 'Van Gogh' travels the longer side.

Your stairlift supplier will advise on which is the best solution for your staircase.

Safety as standard

Your Handicare stairlift will come with safety sensors that are designed to stop your stairlift if it meets any obstruction as you travel up the stairs and will come with a retractable seatbelt as standard.

Easy to use

The operating joystick is easy to use, simply hold it in the direction you wish to travel, up or down, and your Handicare stairlift will move smoothly and gently on its way.